



centri after care

Sales 0800 009 2444

Head Office 01482 329 154

WHW Plastics Ltd.
Therm Road, Hull, HU8 7BF
East Yorkshire, United Kingdom.
Tel. 01482 329154 Fax. 01482 217140
Email. sales@whwplastics.com
Website. www.whwplastics.com
Registered England & Wales 6942043

whw
Behind every smile

NitrAdine

Most removable dental appliances (dentures) are made from porous materials in which opportunistic and pathogenic micro-organisms (fungi, bacteria and viruses) infiltrate, attach and form bio-film, a 3-dimensional protective matrix in which they build-up. Especially, *Candida albicans* penetrates deeply the pores and cracks of the dental appliance. Also other pathological micro-organisms such as *Pseudomonas St mutans*, *St aureus* (inc MRSA), *aeruginosa* and *Herpes simplex* can attach and colonise the appliance which leads to painful irritation of the oral cavity (denture stomatitis).

Elderly, diabetic, transplant and cancer therapy patients wearing contaminated appliances can develop serious systemic infection of the upper airways. For example, oropharyngeal candidiasis is observed in 75% of denture wearing patients in nursing homes. Also the lungs, heart and kidneys can become infected which in some cases may lead to organ failure.

Clinical trials have shown that NitrAdine tablets result in a 99.995% bio-film reduction (Bacteria).

NitrAdine tablets or solution are not to be taken orally

Suitable for flexible dentures

NitrAdine tablet Ingredients (pat pending)

Circic Acid, Sodium Lauryl Sulfate, Lactose Monohydrate, Sodium Bicarbonate, Sodium Chloride, Potassium Hydrogen Monopersulfate, Peppermint Flavour, PVP.

Professional Denture Cleaning Tablets

Our professional effervescent denture cleanser tablets are simple, efficient and easy to use. The combination of the active ingredients removes plaque, impurities, food remains, nicotine, tea, coffee and other stains.

Brushing your dentures alone, is not as effective as a combination of soaking the denture in the effervescent tablets before brushing. By simply brushing, it is easy not to spend enough time cleaning and it is also difficult to clean the hard to reach places which in turn creates nesting places for bacterial growth. By brushing your denture with our NitrAdine Shine Toothpaste and soaking your denture in a special denture solution, you will be assured of giving yourself the best denture hygiene possible.

Are you a **NitrAdine** user? Combine the daily cleaning tabs with a treatment of NitrAdine for the ultimate protection.



Olivafix

A denture fixative cream like no other. What is so special about this product? It's simple. It is the first denture fixative cream on the market that contains **NO Zinc & NO Petrolatum** (more commonly referred to as "Vaseline"). The Zinc and Petrolatum have been replaced by almost 35% Olive Oil. Combining this with a strong, long lasting hold for up to 24 hrs and the fact that Olivafix contains no added preservatives proves that Olivafix is the name for the future for Denture Fixative Creams.

Olivafix is hypoallergenic and leaves a silky feeling in the mouth. It has been clinically tested by consumers with 100% satisfaction and is compatible with all types of dentures. Olive Oil is known for its anti-bacterial and Anti-Candida properties, making it the ideal choice for consumers.

NEW
24h* 'Day & Night' Denture adhesive cream
OLIVAFIX
EXTRA STRONG HOLD
• Long lasting hold
• Neutral taste
• Silky smooth sensation
0% Zinc, 0% Petrochemicals, 0% Parabens

- Strong 24 Hour Hold on full upper denture
- 0% Zinc
- 0% Petrochemicals
- 0% Parabens
- Silky smooth texture
- Neutral taste
- Up to 35% olive oil

Innovative denture
fixative cream

NitrAdine Shine Toothpaste

With its "Anti-biofilm" formula, NitrAdine SHINE has been especially formulated to clean, whiten, refresh and polish your denture or removable appliance. NitrAdine SHINE has been designed to be used on a daily basis as part of your cleansing routine.



Jefe Silicone Denture Brush

Unlike a regular brush, the soft silicone bristles are specifically made for dentures and flexible dentures. The Jefe is non-abrasive and will help keep your denture in excellent condition.



TCS Flexible Denture Cleaner


10 month supply

The perfect companion to your new flexible denture, the TCS flexible cleaner - CLEANS, DISINFECTS & REFRESHES!

Used daily, it will keep appliances looking and feeling like new! Dissolve away stains that brushing leaves behind. Remove coffee, tea, tobacco, and iron stains as well as plaque, tartar, calcium and calculus. It also kills germs that cause bad breath.

Keeps the appliance fresh and clean.

Soak appliance for at least 15 minutes per day. Appliances can also be left soaking overnight or whenever they are not being worn, except appliances containing soft liners, these should not be soaked for more than 15 minutes per cleaning.

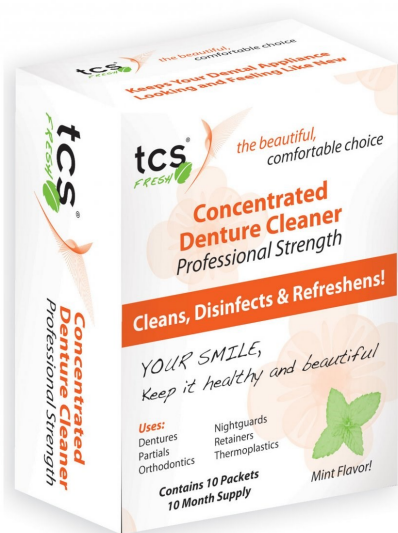

Dent Spa

Denture box with ultrasonic vibrating action for removable appliance cleaning

Insert 2 'AAA' batteries (not included) into the battery compartment located underneath the unit.

Open the unit by pressing the front clip inwards and upwards.

Lightly brush your dental appliance under running water.

Place your dental appliance/s into Dent-Spa and fill with luke warm water just below the lip of the unit (never use boiling water).

Add 1 'NitrAdine' tablet (samples enc.), press the on/off switch to activate. Unit will automatically switch off after 7 minutes.

TABLETS OR SOLUTION NOT TO BE TAKEN ORALLY.

Leave dental appliance/s in Dent-Spa for at least 15 minutes when used with the NitrAdine tablets, then remove appliance/s rinse under water before placing them back in the mouth.

For best results use 'NitrAdine' tablets at least twice a week to leave your appliance fresh, clean and disinfected against harmful bacteria. If you have symptoms of Dental thrush use tablets daily for two weeks to eliminate the bacterial build-up.

Empty Dent-Spa after use, wipe clean, **DO NOT submerge or rinse under running water** as this will damage the unit.

To obtain supplies of 'NitrAdine' or other centri consumables contact either:

1. Your local Dental Surgery
2. Dental Appliance Manufacturer

3. WHW Plastics Ltd
Therm Road
Hull HU8 7BF
Tel 01482 329154

Pk of 20 tablets - NitrAdine(10 wks supply)	£ 10.00	£ .
Pk of 32 tablets - Denture Cleaner	£ 4.75	£ .
1 x 45g NitrAdine Shine Tooth Paste	£ 5.90	£ .
1 x TCS cleaner (40 wks supply)	£ 25.00	£ .
Dent-Spa Case	£ 15.00	£ .
Olivafix	£ 6.75	£ .
Jefe Tooth Brush	£ 5.95	£ .
Post & Packing	£ 1.99	£ 1.99
Please contact your Dental Surgery or Appliance Manufacturer for details on how to make payment		£ .
GOODS SHOULD ARRIVE WITHIN 14 DAYS		

Name.....

Address.....

.....

.....

Post Code.....Tel No.....